

iM3 - The Veterinary Dental Company



iM3 GS Deluxe LED Dental Unit

The iM3 GS Deluxe LED sets the standard for veterinary dental machines.

The GS Deluxe LED is a complete dental system with advanced features not found on any other dental machine.



Quality handpieces coupled with a five-year warranty on the dental machine and oil-free compressor provides peace of mind.

Adjustable stainless steel instrument tray

Advantage Swivel LED High Speed handpiece with excellent ergonomics

Oscillating prophy angles won't catch muzzle hair

Three-way syringe delivering water or CLS, air and mist

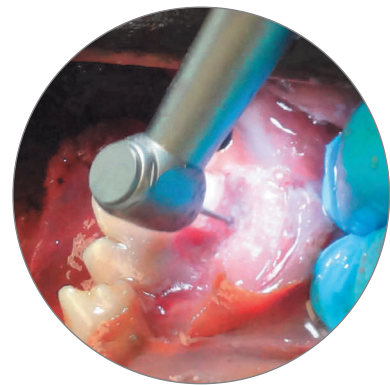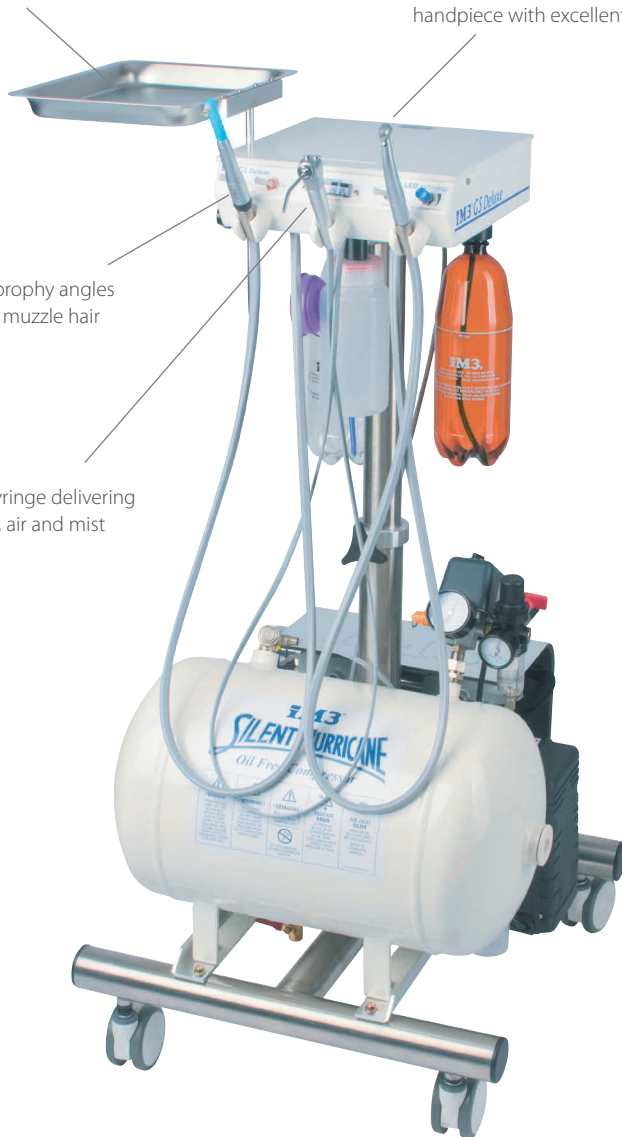

LED High Speed Handpiece Illumination

Light in the mouth is a must!

The new Advantage Swivel LED High Speed Handpiece is the most advanced in the world. It requires no electricity, just air to produce the brightest daylight quality LED white light in the industry with the largest illumination area. Triple-coolant spray provides necessary cooling of the tooth and bur. The patented hygiene head helps to eliminate bacterial suckback and cross contamination.



iM3 GS Deluxe LED with Oil-free Compressor

Cat. No
U3019

Don't forget to add iM3 Straw

Cat. No
PS001

iM3 - The Veterinary Dental Company


Features and Benefits of the GS Deluxe LED Dental Unit



The GS Deluxe LED is a complete dental system with advanced features not found on any other dental machine. The industry's first self-generating Advantage LED Swivel High Speed Handpiece removes the need for troublesome fibre optic wiring, expensive krypton bulbs, switches and transformer. The low speed handpiece and nose cone swivels for comfort and reduces strain on the hand. The GS Deluxe is fitted with two coolant bottles. One for the CLS solution and the other for water. Maintenance instructions are printed on the back of the machine right where you need them.

The GS Deluxe LED is supplied with the following high quality handpieces:

- Advantage Swivel LED High Speed Handpiece
- Three-way syringe (air, water and mist)
- Low speed handpiece

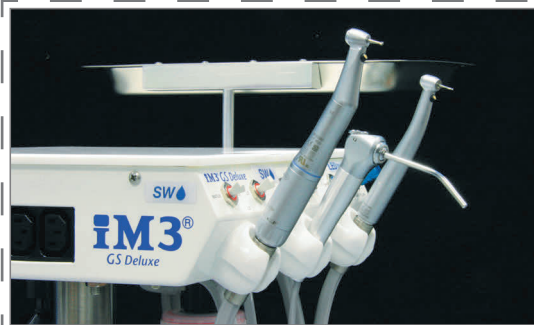


NEW GS Deluxe LED SW (Slow Speed Water) Range

In addition to the standard features found on the GS Deluxe LED, the new SW range includes water on the low speed handpiece.

The SW machine is supplied with an Internal Water Forward & Reverse Low Speed Motor and 1:1 LED Internal Water Low Speed Contra Angle handpiece. Recommended by some European specialists for removing bone.

The SW machine also includes the 4:1 Advantage Straight Nose Cone for polishing.

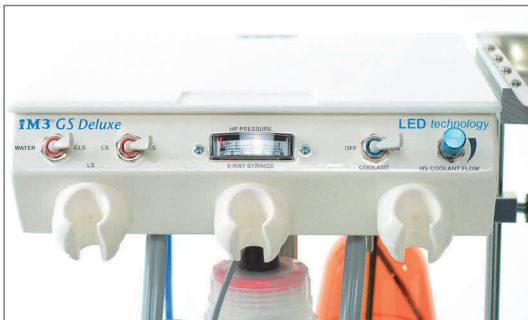

Stainless steel instrument tray provides added storage in the ideal location. It is height adjustable and removable.

The iM3 designed magnetic bur holder prevents burs from falling on the ground; they are held in place right where you need them. The Corian® surface is impervious to bacteria and easy to clean unlike fibreglass and plastics.

The machine features an adjustable coolant flow valve allowing the user to regulate the amount of coolant from the high speed handpiece.

All iM3 toggles, buttons and auto handpiece valves can be changed in seconds due to the cartridge valve system used by iM3. This advanced feature allows for the quick cleaning and replacement of valves reducing machine down time and loss of revenue.



iM3 developed the first oil-free compressor fitted to a dental machine with suction. Oil-free technology eliminates the risk of oil contaminates entering the air supply of the dental machine. It also greatly reduces the maintenance regime. The compressor footprint is 51 x 53 cm (20 x 21"). The oil-free motor generates less heat and is therefore more reliable than conventional oil compressors.

The covered compressor makes cleaning a breeze. iM3's unique technology means the compressor air is cooled on the way to the dental machine. The side of the compressor is left open to allow adequate ventilation, no continuous running fans are required.

Large diameter hospital grade casters make moving the dental machine simple and easy.

The three-way syringe (air, water and mist) is a useful instrument for washing away debris from the mouth. The tip is removable and autoclavable and can be rotated 360 degrees. The use of CLS in the dental machine allows easy application of chlorhexidine into the mouth during dental procedures.



iM3 GS Deluxe LED with Oil-free Compressor

Cat. No
U3019

iM3 GS Deluxe LED SW

(Slow speed with water cooling, not available in all countries)

U3019SW

Don't forget to add iM3 Straw

Cat. No
PS001